

ASMETEC GmbH – Carl-Benz-Str. 4 - D-67292 Kirchheimbolanden – Germany FON: +49-6352-75068-0 – FAX +49-6352-75068-29 – www.asmetec.de - www.asmetec-shop.de – info@asmetec.de

Touptek Microscope - Camera MTR3CMOS USB 3.0 Data Transfer

MTR3CMOS camera adopts SONY Exmor CMOS sensor as the image-picking device and USB3.0 is used as the transfer interface to increase the frame rate.

With the two-stage peltier cooling sensor chip to -42 degree below ambient temperature. This will greatly increase the signal to noise ratio and decrease the image noise.

Smart structure is designed to assure the heat radiation efficiency and avoid the moisture problem.

Electric fan is used to increase the heat radiation speed.

USB3.0 is used as the data transfer interface to increase the frame rate. MTR3CMOS comes with advanced video & image processing application ToupView; Providing Windows/Linux/OSX multiple platform SDK; Native C/C++, C#/VB.NET, DirectShow, Twain Control API;



The MTR3CMOS can be widely used in low light environment and microscope fluorescence image capture and analysis, as well as the astronomy deep sky application.

This camera is available with cylindric housing (standard) or square housing (special order)

The basic characteristic of MTR3CMOS can be summarized as follows:

- Standard C-Mount camera with SONY Exmor CMOS sensors from 1.7M to 45Mpix
- Two-stage TE-cooling with controllable electric fan
- Sensor chip cooling up to 42°C below ambient temperature
- Working temperature can be regulated to specified temperature in 5 minutes
- Smart structure to assure the heat radiation efficiency and avoid the moisture problem
- IR-CUT/AR coated windows
- Up to 1-hour long time exposure
- USB3.0 5Gbit/second interface ensuring high speed data transmission
- Ultra-Fine TM color engine with perfect color reproduction capability
- With advanced video & image processing application ToupView
- Support both video and trigger modes
- Providing Windows/Linux/Mac OS multiple platforms SDK
- Native C/C++, C#/VB.NET, DirectShow, Twain control API;



cylindric (Standard) or square housing

Available versions

Order Code	Sensor & Size(mm)	Pixel(µm)	G Sensitivity Dark Signal	FPS/Resolution	Binning	Exposure
MTR3CMOS45000KMA MTRM145000A	45M/IMX492(M) 4/3 "(19.11x13.00)	2.32 x2.32	702mV with 1/30s 0.12mV with 1/30s	8@8256x5616 31@4128x2808	1x1 2x2	0.1ms~3600s
MTR3CMOS26000KPA MTRP126000A	26M/IMX571(C) 1.8 "(23.48x15.67) APS-C	3.76 x3.76	485mv with 1/30s 0.07mv with 1/30s	6.8@6224x4168(16bit) 14@6224x4168 37@3104x2084 110@2064x1386	1x1 1x1 2x2 3x3	0.1ms~3600s
MTR3CMOS09000KPA MTRP109000A	9M/IMX533(C) 1"(11.28x11.28)	3.76 x3.76	534mv with 1/30s 0.1mv with 1/30s	20@2992x3000(14bit) 40@2992x3000 62@1488x1500 186@992x998	1x1 1x1 2x2 3x3	0.1ms~3600s
MTR3CMOS21000KPA MTRP121000A	21M/IMX269(C) 4/3 "(17.40x13.10)	3.3 x3.3	400mv with 1/30s 0.1mv with 1/30s	5@5280x3956 6@3952x3952 15@2640x1978 50@1760x1318 100@584x440	1x1 1x1 2x2 3x3 9x9	0.1ms~3600s
MTR3CMOS20000KPA MTRP120000A	20M/IMX183(C) 1 "(13.056x8.755)	2.4 x2.4	462mv with 1/30s 0.21mv with 1/30s	5@5440x3648 10@4096x2160 15@2736x1824 30@1824x1216	1x1 1x1 2x2 3x3	0.1ms~3600s
MTR3CMOS20000KMA MTRM120000A	20M/IMX183(M) 1 "(13.056x8.755)	2.4 x2.4	388mv with 1/30s 0.21mv with 1/30s (F8.0)	17.8@5440x3648 41@4096x2160 51@2736x1824 64@1824x1216	1x1 1x1 2x2 3x3	0.1ms~3600s
MTR3CMOS16000KPA MTRP116000A	16M/MN34230PLJ(C) 4/3" (17.6x13.3)	3.8x3.8	2413LSB 89.1LSB (Gain = 0dB)	6@4640x3506 20@2304x1750 48@1536x1160	1x1 2x2 3x3	0.15ms~3600s
MTR3CMOS16000KMA MTRM116000A	16M/MN34230ALJ(M) 4/3" (17.6x13.3)	3.8x3.8	2650LSB 89.1LSB (Gain = 0dB)	22.5@4648x3506 43.0@2304x1750 48.0@1536x1168	1x1 2x2 3x3	0.15ms~3600s
MTR3CMOS10300KPA MTRP110300A	10.3M/IMX294(C) 4/3 "(19.11x13.0)	4.63 x4.63	419mv with 1/30s 0.12mv with 1/30s	7.5@4128x2808 8.5 @4096x2160 30@2048x1080 60@1360x720	1x1 1x1 2x2 3x3	0.15ms~3600s
MTR3CMOS10300KMA MTRM1010300A	10.3M/IMX492(M) 4/3 "(17.47x12.86)	4.63 x4.63	701mv with 1/30s 0.12mv with 1/30s	33.3@4128*2808 69.5@2048*1080 96.2@1360*720 8.0@ 8184*5616	1x1, 2x2, 3x3, Quad*	0.15ms~3600s
MTR3CMOS07100KPA MTRP107100A	7.0M/IMX428(C, GS) 1.1 "(14.4x9.9)	4.5 x4.5	2058mv with 1/30s 0.15mv with 1/30s	12@3200x2200 33@1600x1100	1x1 1x1	0.1ms~3600s
MTR3CMOS07100KMA MTRM107100A	7.0M/IMX428(M, GS) 1.1 "(14.4x9.9)	4.5 x4.5	3354mv with 1/30s 0.15mv with 1/30s	51@3200x2200 133@1600x1100	1x1 1x1	0.1ms~3600s
MTR3CMOS01700KPA MTRP101700A	1.7M/IMX432(C, GS) 1.1 "(14.4x9.9)	9.0 x9.0	4910mv with 1/30s 0.3mv with 1/30s	33@1600x1100	lxl	0.1ms~3600s
MTR3CMOS01700KMA MTRM101700A	1.7M/IMX432(M, GS) 1.1 "(14.4x9.9)	9.0 x9.0	8100mv with 1/30s 0.3mv with 1/30s	94@1600x1100	lxl	0.1ms~3600s
MTR3CMOS08300KPA MTRP108300A	8.3M/IMX485(C) 1/1.2"(11.14x6.26)	2.9x2.9	2188mv with 1/30s 0.15mv with 1/30s	43@3840x2160 66@1920x1080	1x1 2x2	0.1ms~3600s
MTR3CMOS26000KMA MTRM126000A(New)	26M/IMX571(M) 1.8 "(23.48x15.67) APS-C	3.76 x3.76	871mv with 1/30s 0.07mv with 1/30s	6.8@6224x4168(16bit) 14@6224x4168 37@3104x2084 110@2064x1386	1x1 1x1 2x2 3x3	0.1ms~3600s


C:Color; M:Monochrome;

Specifications


Other Specification for MTR3CMOS Cameras				
Spectral Range	380-650nm (with IR-cut Filter)			
White Balance	ROI White Balance/ Manual Temp Tint Adjustment/NA for Monochromatic Sensor			
Color Technique	Ultra-Fine [™] Color Engine/NA for Monochromatic Sensor			
Capture/Control SDK	Windows/Linux/macOS/Android Multiple Platform SDK(Native C/C++, C#/VB.NET, Python, Java, DirectShow, Twain, etc)			
Recording System	Still Picture and Movie			
Cooling System*	Two-stage TE-cooling System -45 °C below Camera Body Temperature			
Operating Environment				
Operating Temperature(in Centidegree)	-10~ 50			
Storage Temperature(in Centidegree)	-20~ 60			
Operating Humidity	30~80%RH			
Storage Humidity	10~60%RH			
Power Supply	DC 5V over PC USB Port External Power Adapter for Cooling System, DC12V, 3A			
Software Environment				
Operating System	System Microsoft® Windows® XP / Vista / 7 / 8 /10 (32 & 64 bit) OSx(Mac OS X) Linux			
	CPU: Equal to Intel Core2 2.8GHz or Higher			
	Memory:2GB or More			
PC Requirements	USB Port:USB3.0 High-speed Port			
	Display:17" or Larger			
	CD-ROM			

Sizes of MTR3CMOS

The MTR3CMOS body, made from tough, alloy with CNC technique, ensures a heavy duty, workhorse solution. The camera is designed with a high quality IR-CUT or AR to block the IR light or protect the camera sensor. The fan's vibration is minimized to the low level to eliminate the vibration caused imaging blur. This design ensures a rugged, robust solution with an increased lifespan when compared to the other industrial camera solutions.

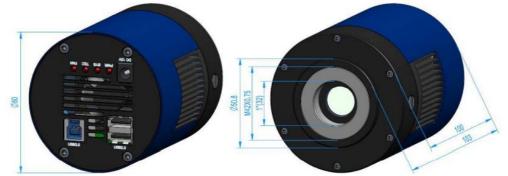

cylindric version (Standard)

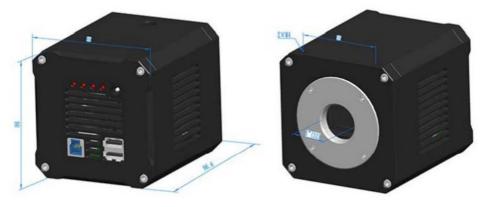

square version (special order)

Packing information



above: cylindric version (standard), below square version (special order)




Supplied with: B, C, D, E, F, G


Adapter H, I, J, K, Glas scalen L are optional

	Standard Package							
Α	Carton L:50cm W:30cm H:30cm (20pcs, 12~17Kg/ carton), not shown in the photo(TBD)							
В	3-A safety equipment case: L:28cm W:23cm H:15cm (1pcs, 2.8Kg/ box); Carton size:L:28.2cm W:25.2cm H:16.7cm(TBD)							
С	MTR3CMOS camera(C-mount)							
D	Drying tube and desiccant							
E	Power adapter: input: AC 100~240V 50Hz/60Hz, output: DC12 V 3A							
F	High-Speed USB3.0 A male to B male gold-plated connectors cable /1.5m							
G	CD (Driver & utilities software, Ø12cm)							
Optional Accessory								
н	Adjustable lens adapter	C-mount to Dia.23.2mm eyepiece tube (Please choose 1 of them for your microscope)	108001/AMA037 108002/AMA050 108003/AMA075 108004/AMA100					
		C-Mount to Dia.31.75mm eyepiece tube (Please choose 1 of them for your telescope)	108008/ATA037 108009/ATA050 108010/ATA075 108011/ATA100					
I	Fixed lens adapter	C-mount to Dia.23.2mm eyepiece tube (Please choose 1 of them for your microscope)	108005/FMA037 108006/FMA050 108007/FMA075 108008/FMA100					
		C-mount to Dia.31.75mm eyepiece tube (Please choose 1 of them for your telescope)	108011/FTA037 108012/FTA050 108013/FTA075 108014/FTA100					
Note: For H and I optional items, please specify your camera type(C-mount, microscope camera or telescope camera), ToupTek engineer will help you to determine the right microscope or telescope camera adapter for your application;								
J	108015(Dia.23.2mm to 30.0mm Ring)/Adapter rings for 30mm eyepiece tube							
K	108016(Dia.23.2mm to 30.5mm Ring)/ Adapter rings for 30.5mm eyepiece tube							
L	Calibration kit	bration kit 106011/TS-M1(X=0.01mm/100Div.); 106012/TS-M2(X,Y=0.01mm/100Div.); 106013/TS-M7(X=0.01mm/100Div., 0.10mm/100Div.)						


Optionale Adapter:

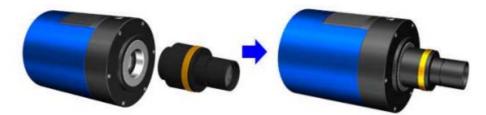


MTR3CMOS+AMAXXX(23.2mm Adapter)

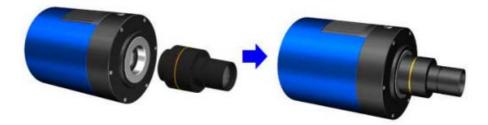

MTR3CMOS+FMAXXX(23.2mm Adapter)

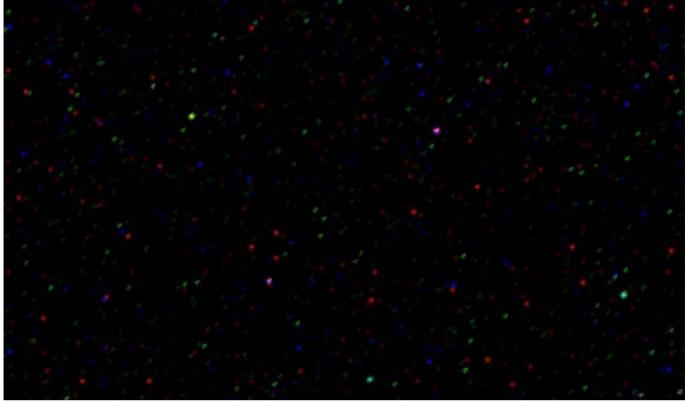

Hot noise for the MTR3CMOS at Gain 20, 600 second, 15 Centidegree

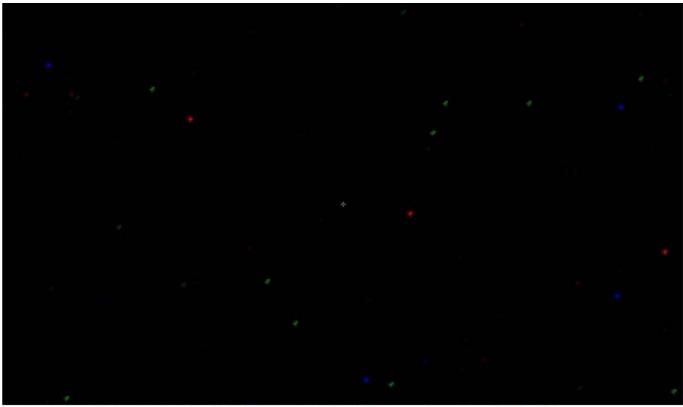

Hot noise for the MTR3CMOS Gain 20, 600 second, minus 15 Centidegree

USB-CAM-MTR3CMOS-DB-E.DOCX, VERSION MRZ-21

The above statements are based on our present knowledge. Our statements should not be interpreted as a guarantee of characteristics. The use of our products by our customers is subject to different conditions, therefore none of our customers are relieved of the responsibility of testing our products by themselves. A liability for consequential damage will not be accepted in any case. For damage resulting from the use of this information we can only be held responsible if there is evidence of malice or negligence on our part. This data-sheet replaces any previous data sheets. ASMETEC, METODRILL, METOCHECK, METOLIGHT, METOSTAT, METOCLEAN and METO are registered trade marks of ASMETEC GmbH.